

Portable Cavitation + Multipolar RF combined machine---SINERON-I

Treatment Scope

- ☑ Soften hard fat tissue, break up lipocyte
- ☑ Strengthen and tighten skin
- ☑ Skin rejuvenation
- ☑ Body shaping, fat reduce, loss weight
- ☑ Cellulite reduce



UNIVERSAL MEDICAL AESTHETICS
Technology designed for your business



Portable Cavitation + Bipolar RF+Tripolar RF +Sixpolar RF combined machine- SINERON-II

Technical Advantages

- ☑ Adopt the top ultrasonic liposuction technique in the world.
- ☑ Suitable for all kinds of skin.
- ☑ Comfortable, painless, woundless during the treatment.
- ☑ Easy operation, easy study.
- ☑ Non- exhaustion, low cost investment.



Products


Portable Cavitation + Multipolar RF combined machine---SINERON-I



Portable Cavitation + Bipolar RF+Tripolar RF+Sixpolar RF combined machine- SINERON-II




Handle function

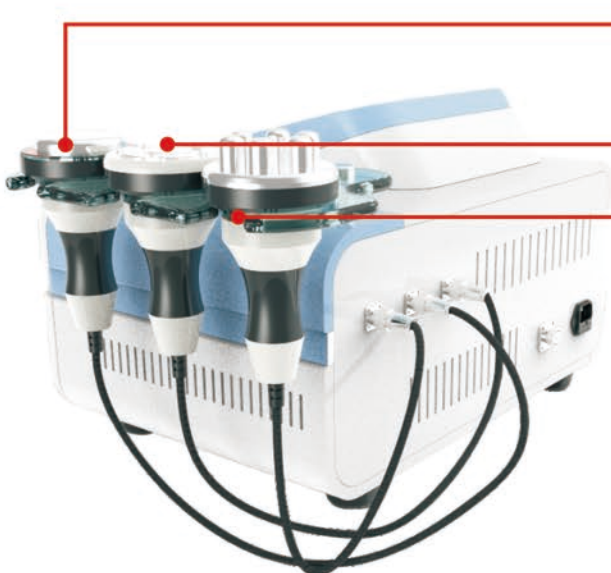

Body sculpting, body shaping, cellulite reduce, etc.

Skin firming, wrinkles improvement, etc.

Skin rejuvenation, Promoting metabolism, Collagen Regeneration

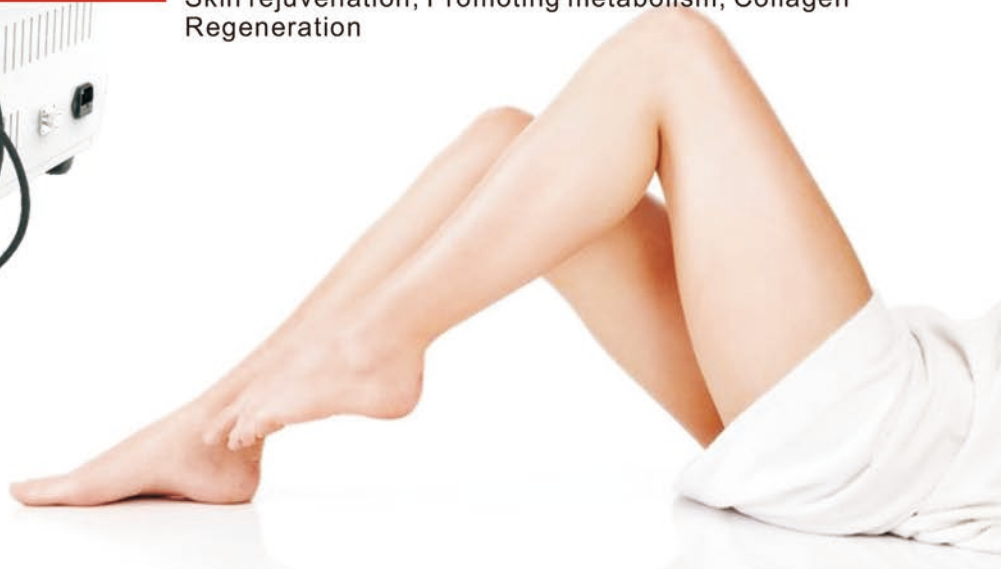

Principle

The Cavitation+RF beauty equipment is new equipment which uses the most advanced technology. It highly combines the advantage of cavitation and RF. The cavitation uses ultrasound cavitation effect to cause the air bubbles around lipocyte membranes to break up. The broken lipocyte will then be absorbed and metabolized out of human body by lymphatic system. It has a perfect treatment result for weight losing. Also we use 5 MHz radio frequency to alternate the electric field in the skin from positive to negative. It will produce heat under dermis then the initial collagen denaturation. Subsequent neocollagenesis then further tightens the dermal tissue and make the skin tightening and smooth.

Specification

Ultrasonic Cavitation Frequency	40KHz	Cavitation Energy	5~30W/cm ²
RF Frequency	5Mhz	RF Energy	0~50J/cm ²
Voltage	110/220V	Power	500W
Package	Aluminum with Foam	Display	10.4 inch digital touch screen
Weight	SINERON-I 20Kg SINERON-II 30Kg	Product size	33cm*31cm*30cm

Before & After

